# ELECTRIC WATER FLOSSER

Item No: 746180

Instructions for Use v001: 14/11/25

Please read and retain these instructions for future reference

## ENGLISH INSTRUCTION

**Power Details:** Input Voltage: 5V

**Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new ELECTRIC WATER FLOSSER.**

### Product Utility


1. The strong pulsed water flow can completely clean up bacteria, tooth stains, tartar and food residues in all parts of the oral cavity.
2. The ultra-fine water stream can clean the gingival sulcus and deep between the teeth without barriers, and improve the health of the gums.
3. It can effectively remove most of the stubborn plaque on the tooth surface, and the effect is much better than traditional dental floss.
4. Quickly clean the oral cavity: Whether at home or eating out, it is the most effective way to quickly restore a clean oral cavity.
5. It has a relatively good effect on the peripheral cleaning of orthodontic teeth, dental implants, and dental veneers
6. Help develop healthy daily interdental cleaning habits and keep teeth fresh and clean at all times.

### Note:

1. Do not use water above 40 degrees Celsius to fill the water tank to flush teeth to avoid burns.
2. Do not use the product for other purposes other than oral cleaning.
3. Do not use the nozzle to press the gums or teeth too hard to avoid damaging the gums.
4. Do not use this product without the nozzle installed.
5. Some functional mouthwashes may be corrosive to a certain extent. Please confirm the ingredients of the mouthwash before use, and do not use it indiscriminately.
6. If children, the elderly or those who do not know how to use this product, please make sure that someone supervises and instructs the user.
7. Please keep it out of the reach of children and babies, and not use it as a toy.
8. Patients with periodontal disease or those who have had oral surgery in the past 2 months, or who have suffered from oral problems, please consult a doctor first and confirm that it is correct before use.
9. Do not use strong gear to continuously wash damaged gums or bleeding teeth to prevent inflammation.
10. Please adjust the flushing intensity of the pulse water flow in real time according to your own oral situation.
11. Do not soak the whole machine in water, it may cause product failure or complete damage.
12. If you do not use this product for a long time, please make sure that the water tank is dry and free of water before storing to prevent the growth of bacteria.
13. The product must be fully charged before use; it must not be used whilst charging.

### Nozzle Specifications:

1. The standard nozzle is suitable for daily cleaning.
2. Tongue scraping nozzle is dedicated to cleaning tongue coating.
3. Periodontal pocket nozzle special for periodontal pocket cleaning.
4. Orthodontic nozzles are suitable for cleaning braces

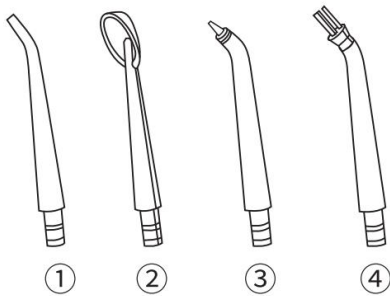
### Method of Operation

1. Long press for 2 seconds to turn on, the default mode is weak, and the corresponding indicator lights up.
2. Short press the switch mode, the mode switches to strong gear, and the indicator light also switches to the corresponding indicator light.
3. Switch the mode to the pulse file, and the indicator light also switches to the corresponding indicator light, and cycle in turn.
4. Press and hold for 2 seconds to shut down. The instrument will automatically shut down after 5 minutes of operation.
5. When the battery of the instrument is charging, the three indicator lights will flash once, and when the battery is full, the three indicators will be on.


### Steps for usage

1. Please fully charge the device before using it for the first time.
2. Insert the corresponding nozzle on the handle to ensure firmness and avoid water leakage during use. (If it be necessary to replace the nozzle, simply pull it out and insert the required nozzle.)
3. Rotate and remove the water tank from the fuselage.
4. Add clean water or mouth wash in the water tank according to your own consumption.
5. Put the nozzle in the mouth and aim at the position to be rinsed.
6. After use, open the back cover to remove residual moisture.
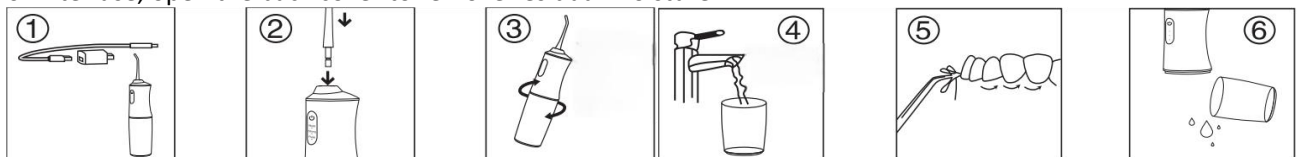
### The right way to clean your teeth:

Choose the gear that suits you. It is recommended to start the trial from the lowest gear and gradually find the gear that suits you best. When cleaning the teeth, extend the nozzle into the mouth, aim at the gap between the teeth or the gums, and open the mouth slightly so that the water flow can flow out smoothly when flushing the teeth. After choosing the gear that suits you, align the gums and rinse them in the same direction from the innermost side (such as 1 or 4), until the gums of the last tooth are rinsed, and then turn the nozzle over, align the gums and continue to the same direction from the inside. Rinse in the direction (such as 2 or 3) until it reaches the gum of the last tooth. The flushing sequence for the other row of teeth is the same as above.
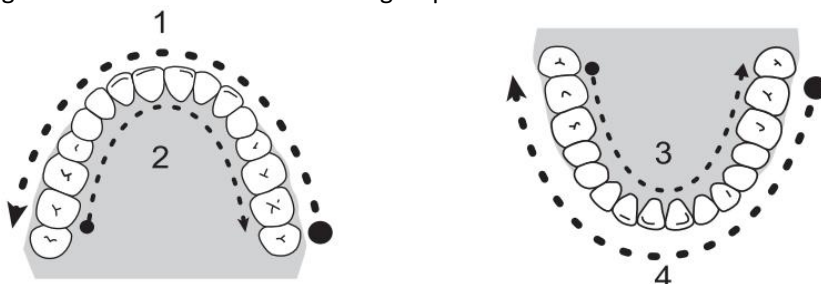
## **WARNINGS AND SAFETY INSTRUCTIONS:**

WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- FOR DOMESTIC INDOOR OR OUTDOORS USE ONLY.
- DO NOT USE FOR COMMERCIAL USE.
- ONLY USE WITH THE USB CABLE SUPPLIED.
- USE ONLY 5V/1A CHARGERS AND CERTIFIED ACCESSORIES. **DO NOT** USE FAST CHARGERS OR INCOMPATIBLE VOLTAGE/CURRENT ADAPTERS AS THIS MAY CAUSE DAMAGE AND INVALIDATES ANY WARRANTY YOU MAY HAVE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- BEFORE CONNECTING THE APPLIANCE, CHECK IF THE VOLTAGE INDICATED ON THE BOTTOM OF THE APPLIANCE IS IN ACCORDANCE WITH THE MAINS VOLTAGE IN YOUR HOME.
- CHECK THE APPLIANCE AND POWER CABLE CAREFULLY FOR DAMAGE BEFORE EACH USE.
- THIS APPLIANCE CAN BE USED BY CHILDREN AGED FROM 10 YEARS AND ABOVE, AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES OR LACK OF EXPERIENCE AND KNOWLEDGE IF THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE IN A SAFE WAY AND UNDERSTAND THE HAZARDS INVOLVED. CHILDREN SHALL NOT PLAY WITH THE APPLIANCE. CLEANING AND USER MAINTENANCE SHALL NOT BE MADE BY CHILDREN WITHOUT SUPERVISION.
- CHILDREN OF LESS THAN 3 YEARS SHOULD BE KEPT AWAY UNLESS CONTINUOUSLY SUPERVISED.
- THE USE OF ANY ACCESSORY OR ATTACHMENTS NOT RECOMMENDED BY COOPERS OF STORTFORD MAY CAUSE INJURIES AND INVALIDATE ANY WARRANTY YOU MAY HAVE.
- DO NOT LET THE POWER CABLE HANG OVER THE EDGE OF A TABLE OR COUNTER, OR TOUCH HOT SURFACES.
- DO NOT PLACE ON OR NEAR HEAT SOURCES.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- ENSURE THERE IS SUFFICIENT DISTANCE BETWEEN THIS APPLIANCE AND OTHER ITEMS ON THE WORK SURFACE OR WORKSTATION, OR SUFFICIENT DISTANCE BETWEEN THE APPLIANCE AND THE CONFINES OF THE AREA IN WHICH IT IS SITUATED. WE RECOMMEND A MINIMUM CLEARANCE AROUND ALL FACES OF AT LEAST 10CM.
- KEEP ITEM AT LEAST 1M AWAY FROM COMBUSTIBLE MATERIAL SUCH AS FURNITURE, CUSHIONS, BEDDING, PAPER, CLOTHES, CURTAINS ETC.
- TO PROTECT AGAINST ELECTRIC SHOCK, DO NOT IMMERSE POWER CABLE, APPLIANCE, IN WATER OR ANY OTHER LIQUID.
- UNPLUG FROM OUTLET WHEN NOT IN USE AND BEFORE CLEANING.
- DO NOT OPERATE ANY APPLIANCE WITH A DAMAGED POWER CABLE OR MAINS ADAPTOR/CHARGER OR AFTER THE APPLIANCE MALFUNCTIONS OR HAS BEEN DAMAGED IN ANY MANNER. WE RECOMMEND THAT THE POWER CABLE IS NOT KINKED OR STORED IN A TWISTED MANNER, AND THAT YOU REGULARLY CHECK ITS CONDITION.
- DO NOT LEAVE UNATTENDED.
- ALWAYS SWITCH OFF AFTER IF YOU HAVE TO LEAVE THE ROOM.

- ANY REPAIRS MUST BE CONDUCTED BY A QUALIFIED ELECTRICIAN OR ELECTRICAL SERVICE AGENT.
- FULLY CHARGE BEFORE FIRST USE.

### CARE AND CLEANING INSTRUCTIONS

- Always switch off and unplug before cleaning.
- Do NOT immerse in water. Risk of electric shock.
- Do NOT use abrasive cleaners or sharp utensils to clean this appliance.
- Wipe surfaces with a soft dry cloth.
- Allow to dry thoroughly before further use.
- Store the appliance in a cool, dry place.

### BATTERIES

- Keep batteries away from children and pets.
- If unused for long periods, recharge every 2 months to protect the battery. Never dispose of batteries in a fire.
- Battery disposal, spent or expired batteries must be properly disposed of and recycled in compliance with local regulations. For detailed information, contact your local authority.
- Follow the battery manufacturer's safety, usage, and disposal instructions.

### DISPOSAL

- We use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons in accordance with your local authority recycling regulations.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.

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- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.

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The Enterprise Department Ltd, Trading as Coopers of Stortford.
**The Enterprise Department Ltd, Trading as – Coopers of Stortford**