VIBRASTONE MASSAGER

Instructions for Use

Item No:692290

v003: 25/03/25

Please read and retain these instructions for future reference

ENGLISH INSTRUCTIONS


Power Details: Power: 12W

Accessories: USB charging cable X 1

Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new vibrastone massager.

INTRODUCTION TO PERATION PANEL:

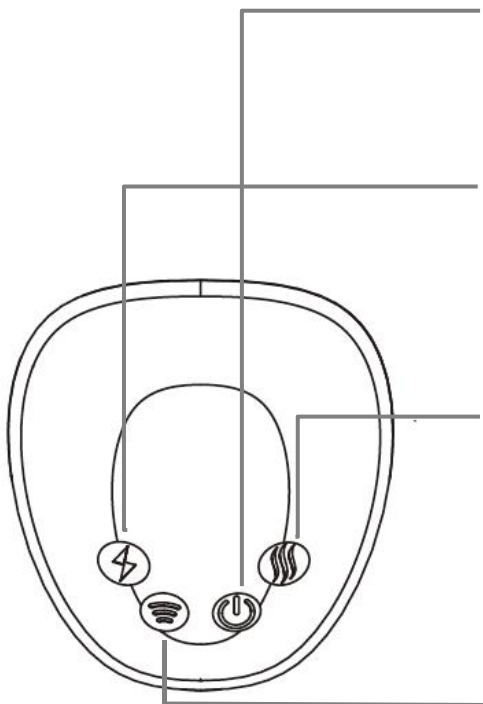Power

Press and hold the button for 3 seconds before turning on the machine, and press and hold the button again for 3 seconds before turning off the machine.

Microcurrent

There are a total of 12 levels of micro current, and after 12 cycles, the micro current function can be turned off. The micro current function needs to be used in conjunction with essential oils for better results. It is recommended to start from a low gear to avoid discomfort caused by excessive micro current flow.

There are a total of 12 levels of thermal energy, which can be turned off after 12 cycles. This function has a heating effect and should be used with caution for people who are not sensitive to heat. The same area should not be used for too long. Please adjust the level according to personal needs to avoid discomfort or burns

Vibration massage

Click on the power button to activate or disable the vibration function. There are a total of 12 vibration levels, and 12 levels are a cycle. After 12 cycles, the vibration function can be turned off

When charging, the indicator light is red and the fully charged indicator light is green. When charging, please use 5V/1A appliances and keep away from flammable and explosive materials such as wood, cotton, paper, etc. Incorrect use of chargers may lead to product malfunctions.

Reminder: The instrument has been set to automatically shut down on a scheduled basis for safety reasons, in order to prevent potential hazards of prolonged heating caused by transportation or accidental operation.

Attention: The product cannot be turned on and used during charging. Please unplug the power to turn it on. When not in use for a long time, please store it fully charged and charge it once a month, otherwise it will damage the battery and make it unusable or shorten its service life.

MASSAGE METHODS

Usage: From top to bottom, massage the left side of the body first, then the right side of the body, starting from the head, neck, back, waist or neck, and back to the limbs;

Head: Use a meridian brush from the front of the forehead to the back of the brain, without the need to apply lubricant randomly. Clear each area about 30 times, and scrape until the head feels hot;

Neck: Brush from top to bottom, you can often scrape and wipe the neck;

Back: Smooth from top to bottom;

Chest: Use a beauty brush to clear the mid line of the chest from top to bottom. The chest is bounded by the Ren meridian on both sides of the anterior mid line of the body;

Abdomen: The abdomen can be cleared from top to bottom (for those with visceral prolapse, it should be cleared from bottom to top), or spiral therapy can be used;

Upper limbs: unblock from top to bottom;

Lower limbs: unclogging from bottom to top (i.e. massaging towards the internal organs);

Sole: Massage from heel to toe.

WARN

This product is not a medical device and is not recommended for medical use.

The product uses electromagnetic pulse (TENS) technology to simulate massage techniques to relax muscles, and any electric flu during the process is normal.

Disabled in the following situations:

1. Patients with skin diseases and skin allergies
2. Pregnant women and individuals in their menstrual period
3. Heart disease patients
4. Patients with abnormal blood pressure
5. Osteoporosis patients
6. Malignant tumor patients
7. Cerebrovascular patients
8. Patients with spinal deformities or abnormalities
9. Individuals with foreign body implants in their bodies
10. Patients with neurological abnormalities in the brain
11. Patients with bleeding tendencies

Please consult a doctor before use in the following situations, otherwise accidents or physical discomfort may occur :

1. Patients with diabetes
2. Patients who cannot express their wishes
3. Patients who are not sensitive to heat
4. Patients receiving medical treatment
5. People who feel unwell
6. Patients with a body temperature exceeding 38 degrees Celsius (fever period)
7. Patients with infectious diseases
8. Patients taking medication
9. Beneficiaries with metal allergies
10. Patients with skin perception disorders
11. Patients with skin abnormalities

TIPS

1. When cleaning the instrument after use, it should not be washed with water. You can use a pure cotton cloth with clean water and wipe it with household alcohol for future use.
2. It is recommended to use each part for about five minutes each time, and a single use time of about ten minutes is recommended.
3. There are 12 adjustable gears for the left micro electric and right thermal energy, and it is generally advisable to turn them to the third to fifth gear.
4. Before using the instrument, apply a hot towel to the area that needs to be massaged, which can moisturize the skin. Because dry and irritable skin can cause uneven penetration, resulting in a stinging sensation.

SAFETY PRECAUTIONS

1. If there is varicose veins in the body, it should not be used.
2. During pregnancy or menstruation, the use should be temporarily suspended.
3. Do not use during charging and keep away from flammable materials.
4. Do not use for a long time, the best massage time is 15 minutes.
5. Children are strictly prohibited from using it, and special users must use it under adult supervision.
6. It is strictly prohibited to dismantle the machine without authorization, otherwise any malfunction or danger will occur, which is not covered by the warranty. Our company will not be responsible for any adverse consequences caused by this.
7. When using, it should be brushed in one direction, not one up and one down, one left and one right.
8. The surface of the appliance generates heat, and people who are not sensitive to heat must pay attention when using it.
9. When massaging, the micro electric brush head should not massage near the heart, otherwise it will increase the risk of heart fibrillation.
10. The product is recommended for one customer to use. When using it for cross users, please clean the brush head with disinfectant before use.
11. This product contains batteries. Please do not disassemble, modify, flood the product, or throw it into fire, as this may cause damage or explosion to the product.
12. If you experience any discomfort during or after use, please stop using immediately and consult a doctor.
13. During massage, the micro electric brush head should avoid contact with wounds and scars.

WARNINGS AND SAFETY INSTRUCTIONS:

WHEN USING ELECTRICAL APPLIANCES, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- FOR DOMESTIC INDOOR.
- NOT FOR COMMERCIAL USE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- THE USE OF ANY ACCESSORY OR ATTACHMENTS NOT RECOMMENDED BY COOPERS OF STORTFORD MAY CAUSE INJURIES AND INVALIDATE ANY WARRANTY YOU MAY HAVE.
- DO NOT PLACE ON OR NEAR HEAT SOURCES.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- DO NOT LEAVE SWITCHED ON OR UNATTENDED.

NOTE

- DO NOT USE VOLTAGE OTHER THAN DC5V FOR CHARGING, AND DO NOT USE IT DURING CHARGING.
- DO NOT USE IN HOT PLACES SUCH AS BATHROOMS, DO NOT IMMERSE YOURSELF IN WATER OR GET CLOSE TO WATER SOURCES TO AVOID SHORT CIRCUITS OR ELECTRIC SHOCK.
- WHEN CLEANING, IN CASE OF MALFUNCTIONS DURING USE, OR DURING USE, THE MACHINE MUST BE TURNED OFF AND THE POWER CORD MUST BE UNPLUGGED.
- IN THUNDERSTORM WEATHER, PLEASE DO NOT CONNECT THE POWER TO CHARGE TO AVOID EXCESSIVE CURRENT, DAMAGE TO THE MASSAGER, AND EVEN PERSONAL INJURY.
- DO NOT CHARGE FOR A LONG TIME. AFTER FULLY CHARGING, UNPLUG THE POWER IN A TIMELY MANNER TO AVOID DAMAGE OR SAFETY ISSUES CAUSED BY OVERHEATING OF THE BATTERY.
- PLEASE USE NEUTRAL CLEANING TO WIPE THE DIAL. DO NOT USE AUTOMOTIVE REPAIR, OIL REFINING, OR HIGHLY ACIDIC OR ALKALINE CORROSIVE CLEANING AGENTS FOR CLEANING.

- AFTER USE, PLACE THE MASSAGER IN A SAFE, DRY, AND COOL PLACE.
- PLEASE AVOID COLLISION OR HEAVY PRESSURE, AND AVOID HARD AND SHARP OBJECTS SCRATCHING THE SURFACE OF THE MACHINE TO AVOID UNNECESSARY MALFUNCTIONS.
- WHEN NOT IN USE FOR A LONG TIME, PLEASE STORE IT FULLY CHARGED AND CHARGE IT ONCE A MONTH, OTHERWISE IT WILL DAMAGE THE BATTERY AND MAKE IT UNUSABLE OR SHORTEN ITS SERVICE LIFE.

BATTERIES


- **BATTERY WARNING: KEEP OUT OF REACH OF CHILDREN. IF SWALLOWED OR PLACED INSIDE ANY PART OF THE BODY SEEK IMMEDIATE MEDICAL ATTENTION.**
- **KEEP BATTERIES AWAY FROM CHILDREN AND PETS.**
- **NEVER DISPOSE OF BATTERIES IN A FIRE.**
- BATTERY DISPOSAL, SPENT OR EXPIRED BATTERIES MUST BE DISPOSED OF AND RECYCLED IN COMPLIANCE WITH YOUR LOCAL AUTHORITY REGULATIONS.
- FOLLOW THE BATTERY MANUFACTURER'S SAFETY, USAGE, AND DISPOSAL INSTRUCTIONS.

BATTERY DISPOSAL

- Never dispose of batteries in a fire.
- Battery disposal, spent or expired batteries must be properly disposed of and recycled in compliance with local regulations. For detailed information, contact your local authority.
- Follow the battery manufacturer's safety, usage, and disposal instructions.

DISPOSAL

- Coopers of Stortford use recyclable or recycled packaging where possible.
- Please dispose of all packaging, paper, cartons in accordance with your local authority recycling regulations.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.
- Plastics, polybags – this is made from the following recyclable plastic.
- At the end of the product's lifespan please check with your local council authorised household waste recycling centre for disposal.

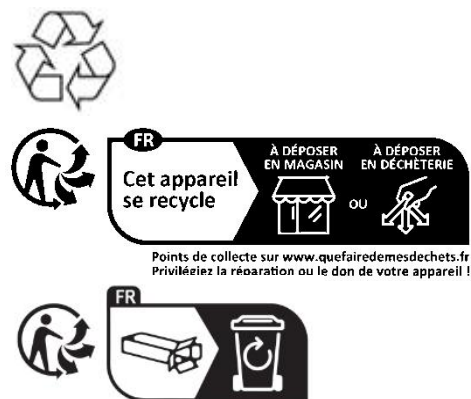Code & Symbol	 LDPE
Type of Plastic	Low Density Polyethylene

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The Enterprise Department Ltd, Trading as Coopers of Stortford.