
GB-Mini Massage Gun

Art. No.: 67117

Massage head

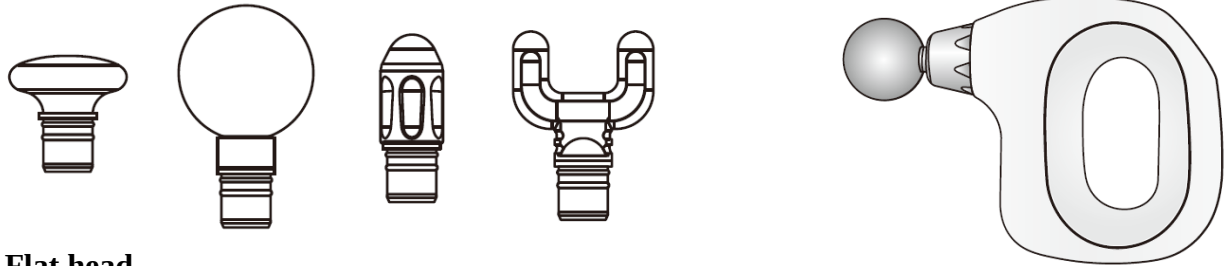-Flat head

Suitable for relaxing and shaping various muscle parts

-Ball head

Suitable for massaging various muscle parts, such as palms and soles

-Cylindrical head

Suitable for impacting deep tissues such as meridian joints

-U-shaped head

Suitable for massaging the neck and spine at the Achilles tendon

Product parameters

Input voltage: 100-240v, 50/60HZ

Operating environment: -5-40°C, non-waterproof

Power supply: 1200 MAH rechargeable lithium battery

Rated power: 20W

PRODUCT USE:

On/off

-Press and hold the start button for 3 seconds to start the machine and enter the first gear.

Force adjustment

-The initial gear for starting the machine is first gear.

-Click the start button once to increase the speed by one gear.

-The maximum gear is twelve, and clicking the start button once at twelve turns it off.

Note: In order to protect bones and muscle groups and avoid long-term fixed point massage that may cause damage, it is advisable to turn on the intelligent protection mode for 15 minutes and automatically stop it.

Indicator light, signal light

- 1、 Tap the switch on and off to display the number, which represents the amount of electricity.
- 2、 Press and hold the key for 3 seconds to start the machine. The default setting for starting the machine is Gear 1. The digital display is 1. The display 2 represents Gear 2, the display 3 represents Gear 3, the display 4 represents Gear 4, the display

5 represents Gear 5, and the display 12 represents Gear 12.

- 3、 When in standby mode, the force becomes the lightest. In energy-saving mode, with slight pressure, the load mode is automatically entered and the force becomes larger.
- 4、 During charging, the digital display number represents the current charging power, and the charging power display is 99.

About Charging

Please be sure to fully charge the battery charging head before using it for the first time.

*To ensure safe use, do not use a fascia gun during charging

*To protect the battery life, please use it fully charged if it is not used for a long time.

Warning:

The following behaviors may cause potential battery safety hazards, leading to safety issues during normal use of the product:

1. Self disassembling battery
2. Firmly press or puncture the product
3. Go to a non-professional maintenance shop for maintenance
4. Use non product standard data lines
5. Place the product in a microwave oven, fire source or near other heat sources.

Security information

1. The product is suitable for working in an environment with a temperature between 0°C and 35°C, and the storage temperature should be between - 20°C and 40°C. High or low temperatures can affect product use and even cause product and battery damage.
2. Please store at temperatures above 5 °C and below 35 °C to avoid degrading battery performance and reducing standby time.

DISPOSAL

-Dispose of all packaging, paper, cartons, plastic and plastic bags in accordance with your local recycling regulations.

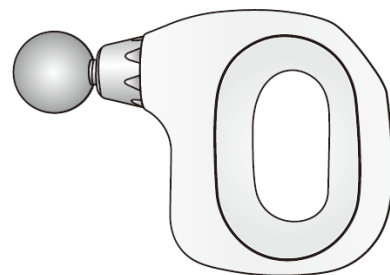
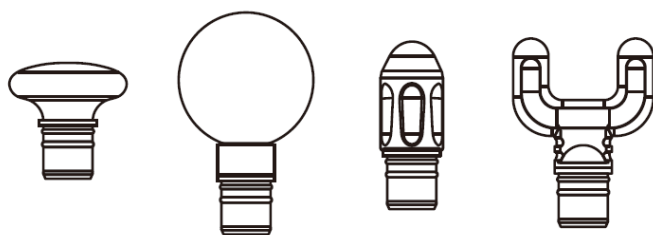
-At the end of the product's lifespan please dispose of it at an authorised household waste recycling centre.

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The Enterprise Department Ltd.