RECHARGABLE ELECTRIC TOOTHBRUSH WITH 4 SPARE HEADS

Item No: 738680

Instructions for Use

v001: 25/08/25

Please read and retain these instructions for future reference

java

wellness


ENGLISH INSTRUCTIONS

Power Details: Battery Type & details: 1200 mAh lithium battery
Voltage details: Input: DC 5V/2A

Accessories: 4 SPARE HEADS + USB charging cable

Please read this instruction manual thoroughly before starting and ensure you are familiar with the operation of your new RECHARGABLE ELECTRIC TOOTHBRUSH WITH 4 SPARE HEADS.

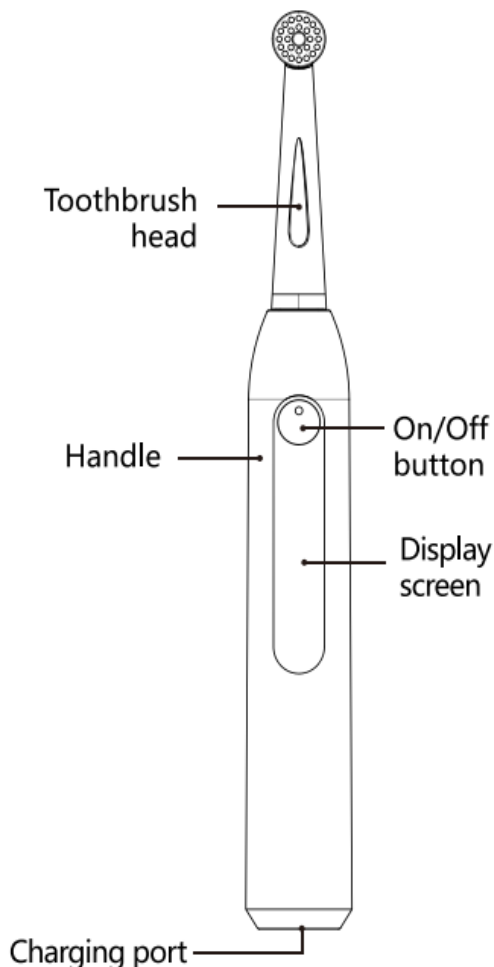CHARGING INSTRUCTIONS: Before the first use, please start turn on the toothbrush. If there is any remaining power, please use it up before charging;

If not used for more than 6 months, please charge the toothbrush for more than 4 hours. The toothbrush will last for 80 days after being fully charged; Please open the bottom cover and insert the charging port of the Type C charging cable into the port of the toothbrush host. The charging indicator "under charge" character on the toothbrush display screen lights up and flashes, and the machine is in charging mode. After being fully charged, the "full charge" character on the toothbrush display screen stays on.

Low battery reminder: The "low power" character on the toothbrush display screen flashes rapidly, reminding you that the battery is about to run out and you need to charge it in a timely manner to avoid affecting use.

If there is no charging for a long time, the indicator light will wait for 6-10 minutes to light up when charging again, which is a normal phenomenon; You should completely discharge it at least once a year, When discharging, unplug the charging cable and let the battery drain normally.

Reminder: This electric toothbrush is suitable for adapters up to 5V-2A and takes around 4 hours to fully charge.

GET STARTED:

1. Insert the brush head to the driveshaft while a certain gap must be left between for the normal vibration.
2. Wet the bristles and apply a certain amount of toothpaste.
3. Put it into your mouth before switching on in case of toothpaste splash. Then press the power bottom again to select a needed mode.
4. Place the bristles against the teeth at 45 degrees towards the gum line.
5. Move the bristles slowly in circular motion. Do not scrub teeth but apply light pressure to let it do the brushing for you.
6. Open the bottom cover and insert the Type C charger plug into the bottom port.

BRUSHING TIME:

The dentist suggests that brushing teeth should last for 2 minutes each time. During use, you will feel a temporary pause every 30 seconds. Here is a reminder that you should clean the four areas of your mouth evenly to achieve the best cleaning effect. If you have slight gum bleeding during the first few uses of your toothbrush, it usually stops within a few days. If the bleeding does not stop, please consult a dentist.

Description of 3 modes:

1. **white:** Loose and stubborn yellow teeth, restore the original white color.
2. **clean:** Strong cleaning power, quickly sweeping away food residues.
3. **soft:** Flexible and gradual, quickly adapting to electric mode.

TROUBLESHOOTING:

If you notice that the vibration has weakened, it may be because the brush head is too tight. Please do not press the bristles too hard.

Reconnect the brush head to the drive shaft, and make sure it is fixed with a certain gap between.

GENERAL NOTICES

This appliance is not intended for use by people (including children) with reduced physical sensory or metal capabilities or lack of experience and knowledge.

Children between the ages of 6-10 must be given supervision or instruction concerning the use of the appliance by a person responsible for their safety. Children under the age of 6 are not recommended for use. This appliance is only for personal use. It was not intended for mixed use in dental clinics by multiple patients.

This appliance is for cleaning teeth, gums and tongue only. Do not use it for other purposes.

Charge indicator may not come to work instantly when you charge after the appliance has been unused for a long period of time. It will switch back on in 6-10 minutes.

It is normal to feel slight heat coming through the handle during the charging process. However, if it turns excessively hot, cut off the power.

Clean the pins of the power adapter after a certain period to avoid overheating and electrocution that may cause by low insulation due to dust accumulation.

Store this appliance away from direct sunlight where it is well-ventilated and dry. Temperature at -7° C ~ 40° C

Do not expose the surface of this appliance to certain substances (gel, oil, soap etc.) for a long time to avoid corrosion.

MEDICAL WARNINGS

It is normal for first-time users to experience slight bleeding and discomfort for their initial week of use. If bleeding or discomfort continues to occur after 1 week of use, discontinue use and contact your dental professional. If you have had oral or gum surgery in the previous 2 months/implanted dental prosthesis, consult your dental professional prior to use.

If you are a pacemaker or other implanted devices, contact your physician to use.

Get related advises if you still have other medical concerns.

WARNINGS AND SAFETY INSTRUCTIONS:

WHEN USING ANY APPLIANCE, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED:

- READ ALL INSTRUCTIONS BEFORE USE.
- FOR DOMESTIC INDOOR USE ONLY.

- NOT FOR COMMERCIAL USE.
- BEFORE USE, CHECK THOROUGHLY FOR ANY DEFECTS AND DO NOT USE IF DEFECTS ARE FOUND. TAKE CARE NOT TO DROP THE APPLIANCE AS HEAVY IMPACTS MAY CAUSE INTERNAL DAMAGE.
- THE USE OF ANY ACCESSORY OR ATTACHMENTS NOT RECOMMENDED BY COOPERS OF STORTFORD MAY CAUSE INJURIES AND INVALIDATE ANY WARRANTY YOU MAY HAVE.
- DO NOT PLACE ON OR NEAR HEAT SOURCES.
- DO NOT COVER THE ITEM IN ANY WAY WHEN IN USE AS THIS MAY CAUSE OVERHEATING.
- DO NOT USE THIS APPLIANCE FOR ANY OTHER USE THAN THE INTENDED USE IT IS DESIGNED FOR.
- DO NOT LEAVE SWITCHED ON OR UNATTENDED.

BATTERIES

- NEVER DISPOSE OF BATTERIES IN A FIRE.
- BATTERY DISPOSAL, SPENT OR EXPIRED BATTERIES MUST BE DISPOSED OF AND RECYCLED IN COMPLIANCE WITH YOUR LOCAL AUTHORITY REGULATIONS.
- FOLLOW THE BATTERY MANUFACTURER'S SAFETY, USAGE, AND DISPOSAL INSTRUCTIONS.

CARE AND CLEANING INSTRUCTIONS

- DO NOT IMMERSE IN WATER.
- DO NOT USE ABRASIVE CLEANERS OR SHARP UTENSILS TO CLEAN THIS APPLIANCE.
- WIPE SURFACES WITH A SOFT DRY CLOTH. DRY THOROUGHLY BEFORE FURTHER USE.
- STORE THE APPLIANCE IN A COOL, DRY PLACE.

BATTERY DISPOSAL

- SPENT OR EXPIRED BATTERIES MUST BE DISPOSED SAFELY AWAY FROM CHILDREN ACCESS AND RECYCLED IN COMPLIANCE WITH YOUR LOCAL AUTHORITY REGULATIONS AND NOT DISPOSED IN GENERAL HOUSEHOLD WASTE.
- FOLLOW THE BATTERY MANUFACTURER'S SAFETY, USAGE, AND DISPOSAL INSTRUCTIONS.

DISPOSAL

We use recyclable or recycled packaging where possible

Please dispose of all packaging in accordance with your local authority's recycling regulations

At the end of the product's lifespan check with your local council authorized household waste recycling centre for disposal.

Plastic Polybags – this is made from the following recyclable plastic

Code & Symbol	 PETE
Type of Plastic	Polyethylene Terephthalate

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The Enterprise Department Ltd, Trading as Coopers of Stortford.